

X-Supreme8000

A powerful, compact, innovative XRF spectrometer delivering highest performance with flexibility



The Business of Science®



Centre of Excellence

For many years Oxford Instruments has been at the centre of innovative science and its practical application worldwide

The first technology business to come out of Oxford University fifty years ago, Oxford Instruments is world renowned for introducing unique, powerful, yet simple to use analytical solutions to a wide range of analysis problems, and this continues with the introduction of the latest X-ray fluorescence spectrometer the **X-Supreme8000**.

Simple, powerful solutions

With the ever increasing requirement for high performance analysis, combined with minimal sample preparation and powerful but simple to use instrumentation, the well respected analytical technique of Energy Dispersive X-ray fluorescence (EDXRF) is finding increasing appeal providing simple analytical solutions to a wide range of quality and process control requirements.

Compared to other analytical techniques EDXRF offers numerous benefits including; non-destructive analysis, minimal or no sample preparation, simultaneous multi-element qualitative and quantitative analysis (results in seconds), and operation by 'anyone'.

Analysis can be performed on a wide range of sample types including solids, liquids, powders, pastes, films, etc, from ppm to high %, covering the elemental range Na¹¹ to U⁹² in the periodic table.

Unique solutions

Oxford Instruments has a long and unique association with EDXRF analysis, for example with the world renowned **Lab-X** and **Twin-X** series of instruments, and to obtain optimal performance involves a close and detailed attention and understanding on the generation and detection of X-rays.

We assure world-class performance in our X-ray tubes and detectors on a global scale. Our Silicon Valley site designs, manufactures and distributes high-reliability X-ray tubes for both industrial and medical applications. We are also heavily involved in the European space program with the design and production of solid-state X-ray detectors. These technologically superior detectors, developed in our Finland production facility, are used extensively in our industrial XRF analysers.



Image: Oxford spires

The Business of Science[®]

Introducing the X-Supreme8000

The **X-Supreme8000** is a compact, powerful, yet easy to use EDXRF spectrometer. It incorporates an application optimised Oxford Instruments' high reliability X-ray tube and high-performance Silicon Drift (SDD), or proportional detector.

This unique best-in-class combination gives highest performance with lowest detection limits, and provides a winning integrated solution for quality assurance and process control requirements.



Cement



Wood treatment



Minerals



Mining



Cosmetics



General analysis



Petrochemicals



Education

Supreme technology for high performance and flexibility in XRF analysis



X-Supreme8000 As easy as...

1

Prepare sample



Minimal or no sample preparation

- No weighing or volumetric measurements
- No potentially dangerous chemicals
- Non-destructive analysis

2

Place sample on tray



Easy loading

- 10 positions for time saving unattended operation
- Unique design guarantees high precision, consistent, repeatable results from every tray position
- Tested to 1 million cycles, fully loaded – equivalent to 15 years of operation at 200 samples per day

3

Press start key



Rapid results

- Simultaneous single or multi-element analysis with 'Live Updates' – results obtained 5 seconds after measurement starts, allowing a rapid assessment to be made
- Results obtained using empirical calibrations for highest accuracy and traceability. Optional Standardless (Fundamental Parameter) Analysis where no or very few standards are available
- Rugged and reliable operation on a 24/7 basis carried out by laboratory or shift personnel



**Results – save,
print or send!**



Speed, accuracy and long term reliability – FocusSD technology

Highlighted area shows excitation and detection – FocusSD technology in action



- Application optimised, high performance solutions
- Oxford Instruments' X-ray tubes with W, Pd or Ti target and state of the art high resolution large area Silicon Drift (SDD) or proportional detector

Minimum bench space

- Integrated PC – compact
- Large, easy to read colour graphics for clarity and ease of use
- Touch screen display (optional) for simple and rapid data input
- High performance, reliable, embedded PC with Windows XP

Easy communication

- Easy access USB ports on side/underneath PC module for peripheral devices
- Ethernet connection for; transmission of results, interface to Laboratory Information Management System (LIMS) or remote operation and diagnostics



Unattended operation



- 10 sample positions with ability to measure QC sample during run giving total confidence in analytical results
- Sample spinner (optional) – reduces errors due to sample preparation, particle size and inhomogeneity
- Optional large trays for holding 51.5 mm (2") diameter steel ring for cement application, or 47 mm diameter filter samples
- Range of Oxford Instruments' sample holders for easy sample loading

Rugged and robust design

- Sealed membrane keypad resistant to dust, oils, solvents etc.
- Internal 'wind tunnel' design – cooling fans do not carry dust to sensitive internal components
- Secondary easy-to-change safety window for liquids/powder analysis protects the instrument in the event of sample leakage


Comprehensive analysis

The **X-Supreme** software uniquely provides simplicity of operation for routine analysis, with flexibility for qualitative, semi-quantitative and full quantitative analysis.

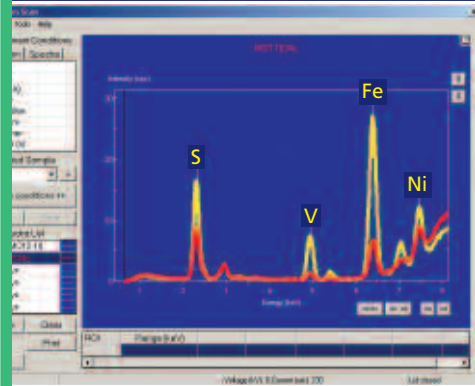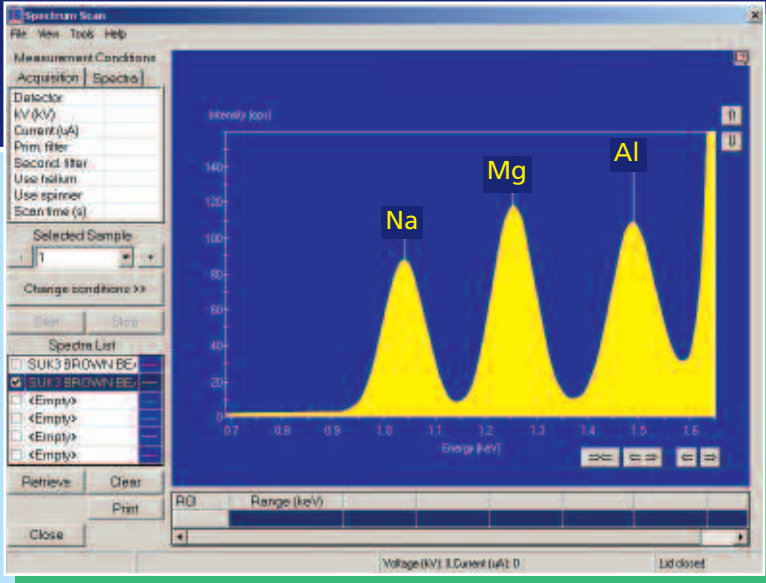


Investigation (qualitative analysis)

- Focus^{SD} – unique elemental sensitivity and separation
- Multi-element X-ray spectrum obtained in seconds
- High resolution, high performance

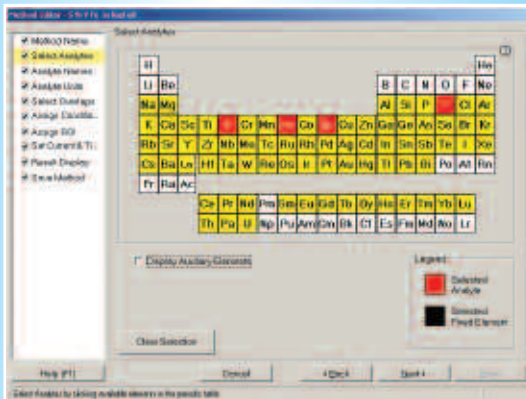

X-Supreme method Wizard

- Simple step-by-step setup
- Straight forward editing

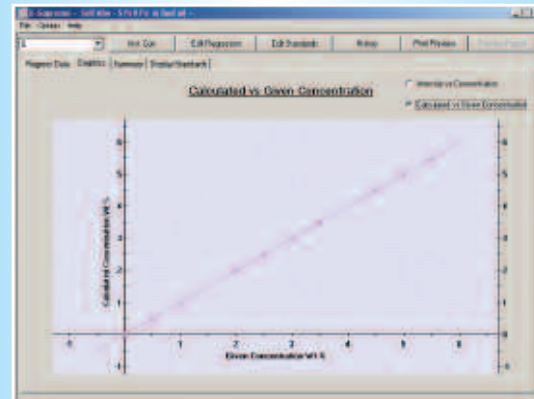


Calibration (Quantitative Analysis)

- Regression – intensity, ratio, background corrections etc.
- Highest accuracy and traceability of analysis



Routine Analysis (for production staff)

- Helpful hints for advice on good practice
- QC sample measurement ensuring high quality analysis
- SMARTCHECK software allowing operation by anyone

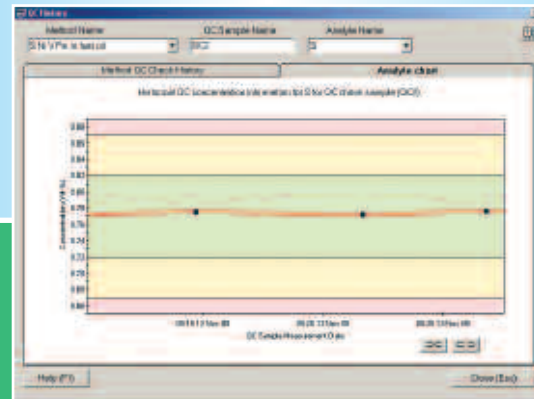

Quality control

- Monitoring of instrument performance
- Powerful graphical/numerical QC display ensures high analytical performance

Worldwide customer service and support

Oxford Instruments takes great pride in partnering with our customers to provide comprehensive technical support and instrument service throughout the product's life.

The Oxford Instruments Customer Support Team delivers the background technical maintenance, expertise and service to both customers and our global distributor network to keep your analyser running smoothly.

Our global network provides a wide range of customer support functions:

- Basic and Advanced user courses
- Support and service contracts for minimum downtime, and ISO audit and certification purposes
- Oxford Instruments' expert technical helpdesk
- Warranty support and instrument repair
- Instrument (software and hardware) upgrades
- Spare parts (available through trained distributors)
- Consumables – complete alphabet of products ranging from Aluminium sample cups to Zirconium standards



Behind your success

Your investment in our equipment is backed by global and local support of the highest standards, with the aim of providing a comprehensive support package to ensure your success.

click onto www.oxford-instruments.com for more information

Oxford Instruments, at High Wycombe, UK, operates Quality Management Systems approved to the requirements of BS EN ISO 9001. This publication is the copyright of Oxford Instruments Analytical Limited and provides outline information only which (unless agreed by the company in writing) may not be used, applied or reproduced for any purpose or form part of any order or contract or be regarded as a representation relating to the products or services concerned. Oxford Instruments' policy is one of continued improvement. The company reserves the right to alter, without notice, the specification, design or conditions of supply of any product or service. Oxford Instruments acknowledges all trade marks and registrations.

© Oxford Instruments Analytical Ltd, 2009. All rights reserved.

As part of Oxford Instruments' environmental policy this brochure has been printed on FSC paper.



Certificate No FM29142

Part no: OIIA/069/A/0109

Oxford Instruments Industrial Analysis

UK

Halifax Road, High Wycombe
Bucks, HP12 3SE England
Tel: +44 (0) 1494 442255
Fax: +44 (0) 1494 461033
Email: analytical@oxinst.com

China

Beijing
Tel: (8610) 6518 8160/1/2
Fax: (8610) 6518 8155
Email: info@oichina.com

Finland

Espoo
Tel: +358 9 329 411
Fax: +358 9 3294 1300
Email: FI-Espoo_Info@oxinst.com

Germany

Uedem
Tel: +49 (0) 2825 93 83 -0
Fax: +49 (0) 2825 93 83 -100
Email: DE-Uedem_Info@oxinst.com

Japan

Tokyo
Tel: +81 (0) 3 5245 3591
Fax: +81 (0) 3 5245 4466/4477
Email: oikkma@oxinst.co.jp

Latin America

Clearwater FL
Tel: +1 727 538 7702
Fax +1 727 538 4205
Email: oxford@gate.net

Singapore

Tel: +65 6337 6848
Fax: +65 6337 6286
Email: asiasales@oxinst.com

North America

Elk Grove Village IL
Tel: +1 847 439 4404
Fax: +1 847 439 4425
Email: sales@msys.oxinst.com

www.oxford-instruments.com

